



User Guide : Epoxym Kit

Create Robust Replica Molds

EPOXYM 4

DESCRIPTION	5
EPOXYM COMPONENTS	5
EQUIPMENT	5
SAFETY NOTES	6

PROCESS 7


OPERATING INSTRUCTIONS	8
EPOXY CALCULATION REFERENCE TABLE	13
MACHINE MAINTENANCE	13

TERMS & CONDITIONS 14

DESCRIPTION

EpoxyM is a user-friendly solution designed for replicating SU8 master wafers and other types of molds. The replication process involves two steps: first, creating a polydimethylsiloxane (PDMS) positive replica, and then using this as a template to produce an epoxy resin replica of the original master.

EPOXYM COMPONENTS

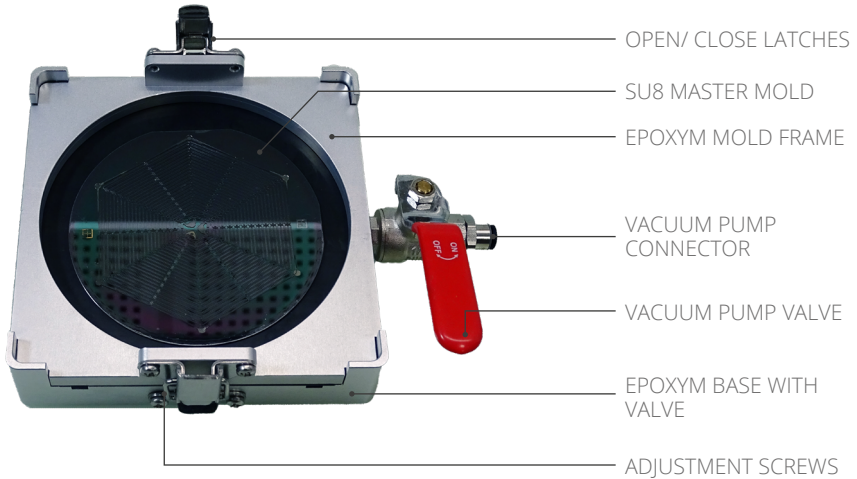

EQUIPMENT INCLUDED IN THE KIT

- EpoxyM base with valve
- EpoxyM mold frame for PDMS (Ø: 12 cm)
- EpoxyM mold frame for Epoxy (Ø: 10 cm)
- O-RING, FKM, 126 x 3 mm
- Flexdym™ seal sheet
- Epoxy (resin & hardener)
- Silicon spray
- Open/close latches with adjustment screws x2

EQUIPMENT NOT INCLUDED IN THE KIT

- PDMS (base & curing agent)
- Master mold
- Desiccator
- Vacuum pump with 6 mm OD tubing
- Ethanol or isopropanol
- Master mold (Ø: 4")
- Vented oven or Hot plate

SAFETY NOTES

The aim of this quick start guide is to illustrate how to use the Epoxym equipment for Epoxy mold replica fabrication.

Please, always consult the safety notes or procedures that may be related to your specific equipment prior to use.

As with all laboratory activities, gloves and lab coats should be worn at all time during the process. It is the user's responsibility to consult the manufacturer's recommendation for the chemicals used and to ensure adequate ventilation/extraction (e.g., via a fume hood) prior to proceeding. Do not exceed the maximum temperature stated in the equipment specification.

Epoxy handling:

Always use a mask, gloves and protective clothing.

Avoid breathing mist or vapours.

Wash skin thoroughly in case of contact.

Avoid release in the environment.

Epoxy resin (Part A):

Causes skin irritation.

May cause an allergic skin reaction.

Causes serious eye irritation.



Hardener (Part B):

May cause an allergic skin reaction.

Causes serious eye damage.

May cause allergy or asthma symptoms or breathing difficulties.



OPERATING INSTRUCTIONS


The Epoxym kit enables the production of durable epoxy molds resistant to high temperatures, up to 180°C. The process begins with the conversion of an SU-8 master mold into an intermediate PDMS positive replica, which is then used to produce the epoxy mold.

1. Preparation of PDMS

Prepare the PDMS solution by mixing the base & hardener according to the manufacturer's instructions (e.g. ratio of 1:10 for Sylgard™ 184). Evaluate the PDMS volume to get 4-6 mm thick PDMS positive replica.


Desiccate the PDMS mix for at least 30 minutes to remove all air bubbles.

Tip: For Sylgard™ 184, 60g of base and 6g of curing agent should be sufficient.

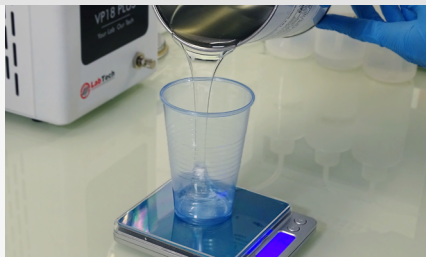

2. Preparation of the Master Mold

In the meantime, gently clean the silicon master mold to be replicated with isopropanol and preferably use an air gun to dry it.

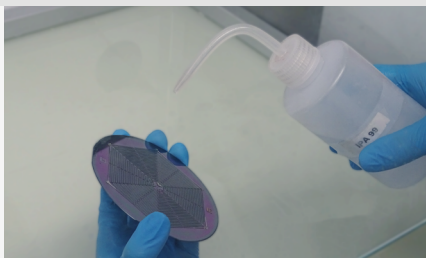


3. Preparation of the Epoxym Kit

Place the Flexdym™ black sealing sheet on the base and center the master mold on top.

Insert the rubber O-ring seal into the PDMS frame, then fix it to the base closing the side latches. Connect the setup to the vacuum pump, switch it on, and open the valve.


Check that the mold is securely fixed before proceeding.

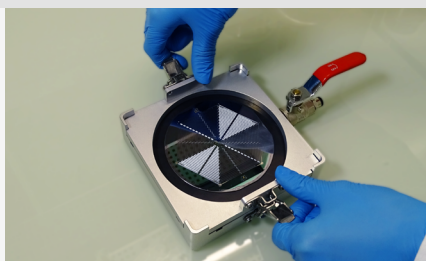

4. Pour PDMS

Very gently, pour the degassed PDMS into the mold frame.

Note: This process should be done slowly to not produce air bubbles at the level of the microstructures.

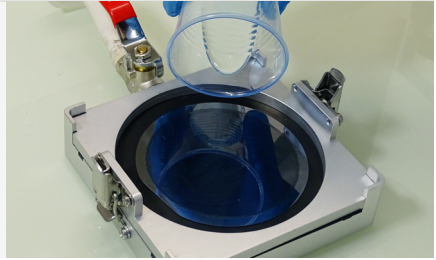


5. Bake PDMS

Close the vacuum valve and disconnect the Epoxy from the pump.

Bake the PDMS in the oven according to manufacturer's instructions (e.g. 2 hours at 80°C for Sylgard™ 184).

Note: Make sure the oven shelf is flat.

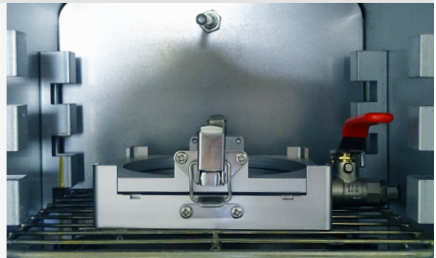


6. Unmold PDMS replica

Remove the Epoxy from the oven and unmold the PDMS from the frame.

Open the vacuum valve and remove the master mold and Flexdym™ sealing sheet from the base using flat tweezers.

Tip: Protect the PDMS microstructures from dust by covering them with adhesive tape.

Note: Store the Flexdym™ sealing sheet in a clean environment to protect it from dust and maintain its adhesive properties.

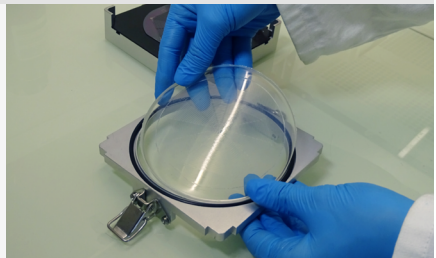

7. Preparation of Epoxy Resin (1/2)

Weigh the epoxy components A (resin) and B (hardener) according to the calculation table (page 14).

Heat up the components separately on a hot plate at 70°C for at least 10 minutes.

Warning: Always perform epoxy handling and heating in a fume hood or vented oven.

Note: After use, clean the edges of the bottles thoroughly to prevent dried epoxy from sealing them shut. Seal the bottle caps with parafilm and store at room temperature.




8. Preparation of Epoxy Resin (2/2)

Mix components A and B of the epoxy using a stirrer.

Heat at 70°C for another 10 minutes.

Tip: By pouring component B (hardener) into component A (resin), volume waste is minimized.

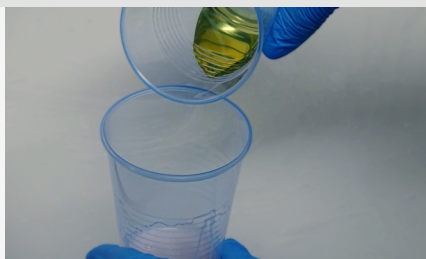


9. Preparation of the Mold Frame

Coat the epoxy mold frame with silicone spray to allow the unmolding step.

The side latches should be screwed into the threaded slot that corresponds to the thickness of the PDMS counter-mold.

Adjust the latches position to ensure proper tightness and an effective seal of the frame.

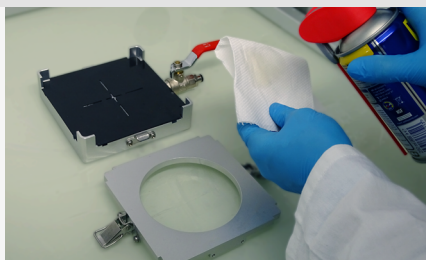

10. Preparation of the Epoxy Kit


Place PDMS replica face up in the center of the base and use the epoxy mold frame to close the system via the side latches.

Pre-heat the Epoxy system at 70°C, then connect to vacuum pump.

Open valve and apply vacuum. Remove the tape from the PDMS replica.

Tip: Pour a small amount of isopropanol onto the PDMS surface and rotate the Epoxy system to verify the tightness and seal of the frame.


Note: The PDMS replica must be flat. If a curvature is observed, it may result from a "meniscus" formed at the edges. Carefully remove this excess material with a scalpel.

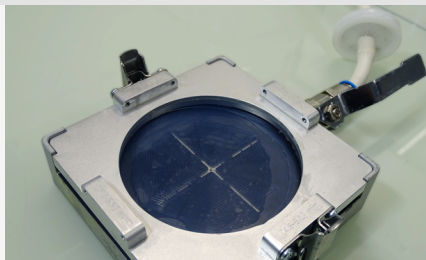

11. Cast Epoxy on PDMS replica

Pour the epoxy mix into the mold frame very gently, in order to not create bubbles.

Note: If the structures to be replicated are small and intricate, place the Epoxy system (with the vacuum valve open) in the dessicator to remove eventual air bubbles from the resin.

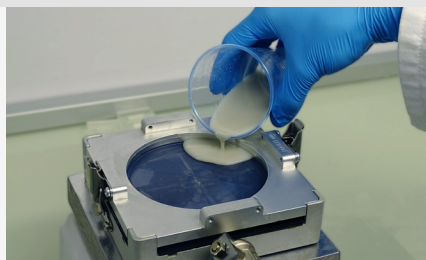

12. First Bake

Place the Epoxy - still connected to the vacuum pump - on a hot plate under the fume hood. Make sure it is flat and bake at 90°C for 3 hours, 120°C for 2 hours, and at 150°C for 2 hours.


Tip: Keep the vacuum pump connected to the Epoxy base until the resin fully solidifies. This prevents buckling of the PDMS replica, which could compromise the flatness of the final epoxy mold.

Note: This step can be performed in a vented oven. If using an oven, first remove the adapter that connects the Epoxy base to the vacuum pump tubing, as it contains plastic components that may melt. Ensure the oven shelf is flat before proceeding.




13. Unmold

Open the side latches and unmold the epoxy mold from the frame.
Cover the PDMS counter-mold with a scotch tape for later use.

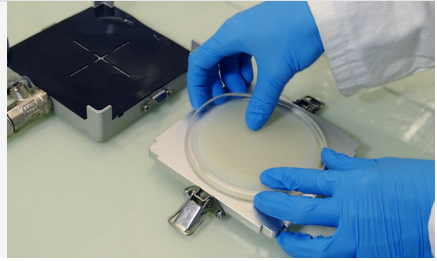

14. Hard Bake

Hard bake the epoxy mold at 170°C on the structure side and 180°C on the bottom side for 1 hour.
To ensure the mold stays flat, we suggest using a heating press or to use a hotplate with a weight on top of the mold.

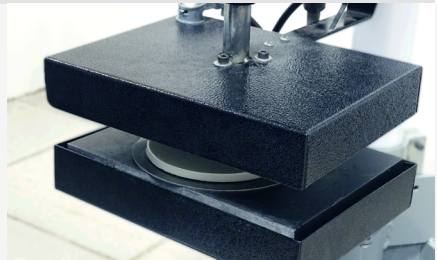

EPOXY CALCULATION REFERENCE TABLE

The recommended epoxy mold thickness is 3 to 6 mm.

To achieve this, approximately 35 ml of resin should be sufficient.

The table below indicates the quantities to be used for each epoxy resin component:

	COMPONENT A : RESIN	COMPONENT B : HARDENER
PARTS BY VOLUME	100	67
PARTS BY WEIGHT	100	83

MACHINE MAINTENANCE

After each use, we recommend cleaning the aluminum parts with acetone and ethanol, and the O-ring with ethanol only. Ensure the frame is coated with silicone spray before each use, and promptly clean any epoxy spills before polymerization to maximize the shelf life.

Avoid contact between metallic tools and the frames, as it might damage their surface finish and promote leaks.

TERMS & CONDITIONS

We strongly believe in the intrinsic quality of our microfluidic products and we hope that you will be pleased with your purchase. However, in the unlikely event that you should receive damaged or incorrect goods, please notify us within 7 days of receipt and you will be offered the option of a refund or an exchange (provided the goods are in stock). You may be asked to return goods for inspection. In this case we will refund the shipping fees. Please note that goods that become damaged or broken after 7 days of receipt cannot be returned.

Unwanted items

If for any reason you do not wish to keep your purchase and would like a store credit, then please notify us within 7 days of receipt. We cannot accept unwanted returns that have been opened, used or damaged by the customer. For unwanted goods, we allow up to 14 days for the returns of goods. We will only issue a credit upon receipt of all returned goods.

Privacy policy

Customer details remain private and confidential and will not be released to a third party unless required to do so by law. We use the information we collect to process orders, to provide more personalized shopping experience and, if you request it, to notify you about new products, special offers or other information that may be of interest to you.

Payment and Procedure

Payment from private customers must be paid for in advance of shipment. Trade orders from registered companies or organization can be invoiced. Payment is due strictly within 30 days of the invoice date.

Products and prices

Please note that some goods may vary in style, color or detail from the image shown online. We reserve the right to change prices at any time.

Transport and storage

Please prevent damage to your Eden Tech products by taking extra care when moving. Products must not be transported while plugged in or hot. Store products in standard room condition and in an adapted box.

Exclusive remedies

The remedies provided herein are the customer's sole and exclusive remedies. Eden Tech shall not be liable for any direct, indirect, special, incidental, or consequential damages, whether based on contract, tort, or any other legal theory.

Warranty

The EpoxyM Kit is covered by a one-year warranty from the date of purchase. This warranty ensures that the machine is free from defects in materials and workmanship under normal use and conditions. Any issues arising during the warranty period will be assessed, and repairs or replacements will be provided as necessary. The warranty does not cover damage caused by improper use, unauthorized modifications, or external factors beyond our control.



Eden Materials is a division of Eden Tech

Registered Address: 83 Avenue Philippe Auguste, 75011 Paris, France | Client Facing Address: 4 rue de Rambervillers, 75012 Paris, France

T. + (33) 188 327 405 E. contact@eden-microfluidics.com W. www.eden-microfluidics.com